



PB Hemo2Way

INDICATIONS

The PB Hemo2Way is designed to perform coagulation, dissection, cutting, or electrosurgical fulguration, while providing continuous fluid irrigation and suction through the Viewing Cannula in delicate surgical procedures. The irrigation is done by gravity, but the flow can be amplified with the use of the unique Flow Booster, as well as easily controlled by the surgeon with the clamp to interrupt the flow. Its active tip in noble and polished material avoids tissue adhesion and provides heat transfer away from the tip, avoiding overheating of adjacent tissues.

Its active tip, which reproduces the natural pinching movement of the hands, facilitates hemostasis, offering safety and precision in the surgical act. Suction and continuous irrigation favor a clean field of vision that results in little bleeding in highly vascularized regions. Such characteristics allow the preservation of the structures, avoiding possible bleeding and improving the post-operative period.

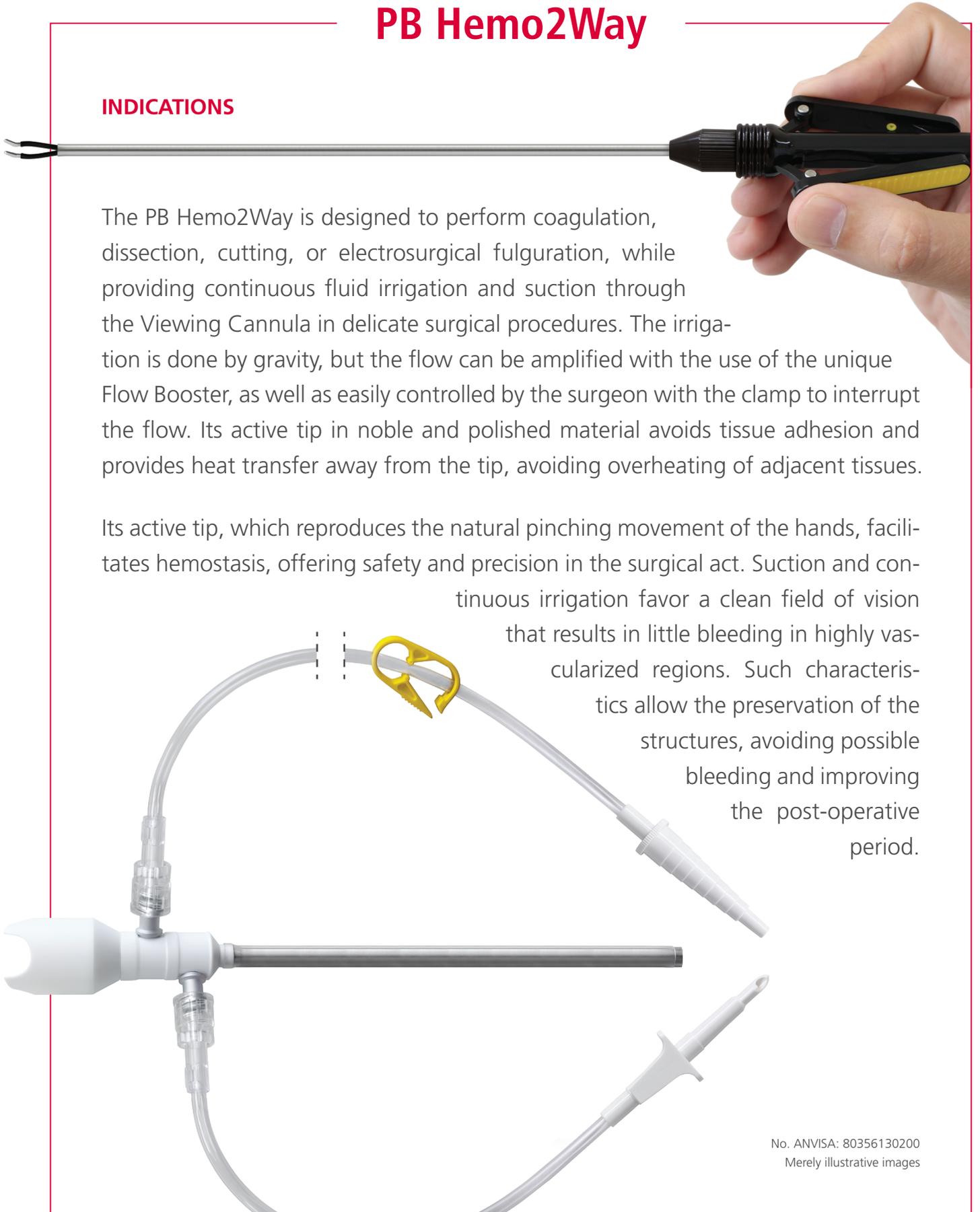

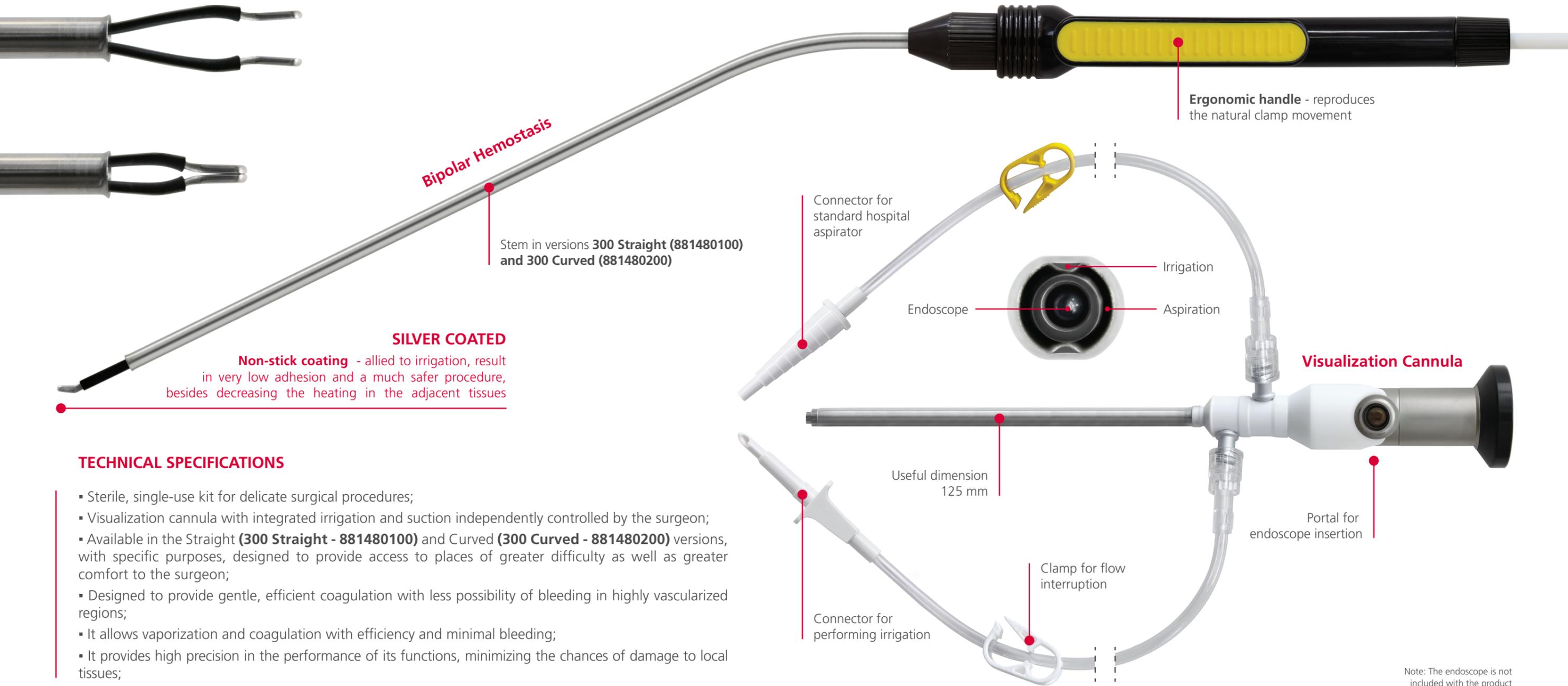

TECHNICAL SPECIFICATIONS

- Sterile, single-use kit for delicate surgical procedures;
- Visualization cannula with integrated irrigation and suction independently controlled by the surgeon;
- Available in the Straight (**300 Straight - 881480100**) and Curved (**300 Curved - 881480200**) versions, with specific purposes, designed to provide access to places of greater difficulty as well as greater comfort to the surgeon;
- Designed to provide gentle, efficient coagulation with less possibility of bleeding in highly vascularized regions;
- It allows vaporization and coagulation with efficiency and minimal bleeding;
- It provides high precision in the performance of its functions, minimizing the chances of damage to local tissues;

Note: The endoscope is not included with the product



- Bipolar hemostasis clamp with an active tip that allows for tissue clamping and an ergonomic handle that reproduces the natural clamping movement of the surgeon's hand.

ACCESSORY



Flow Booster: The accessory must be used on the fluid container to intensify the flow of irrigation fluid. 2 units per product are included.





www.razek.com.br
+55 16 2107.2345